


Polypectomy Solutions



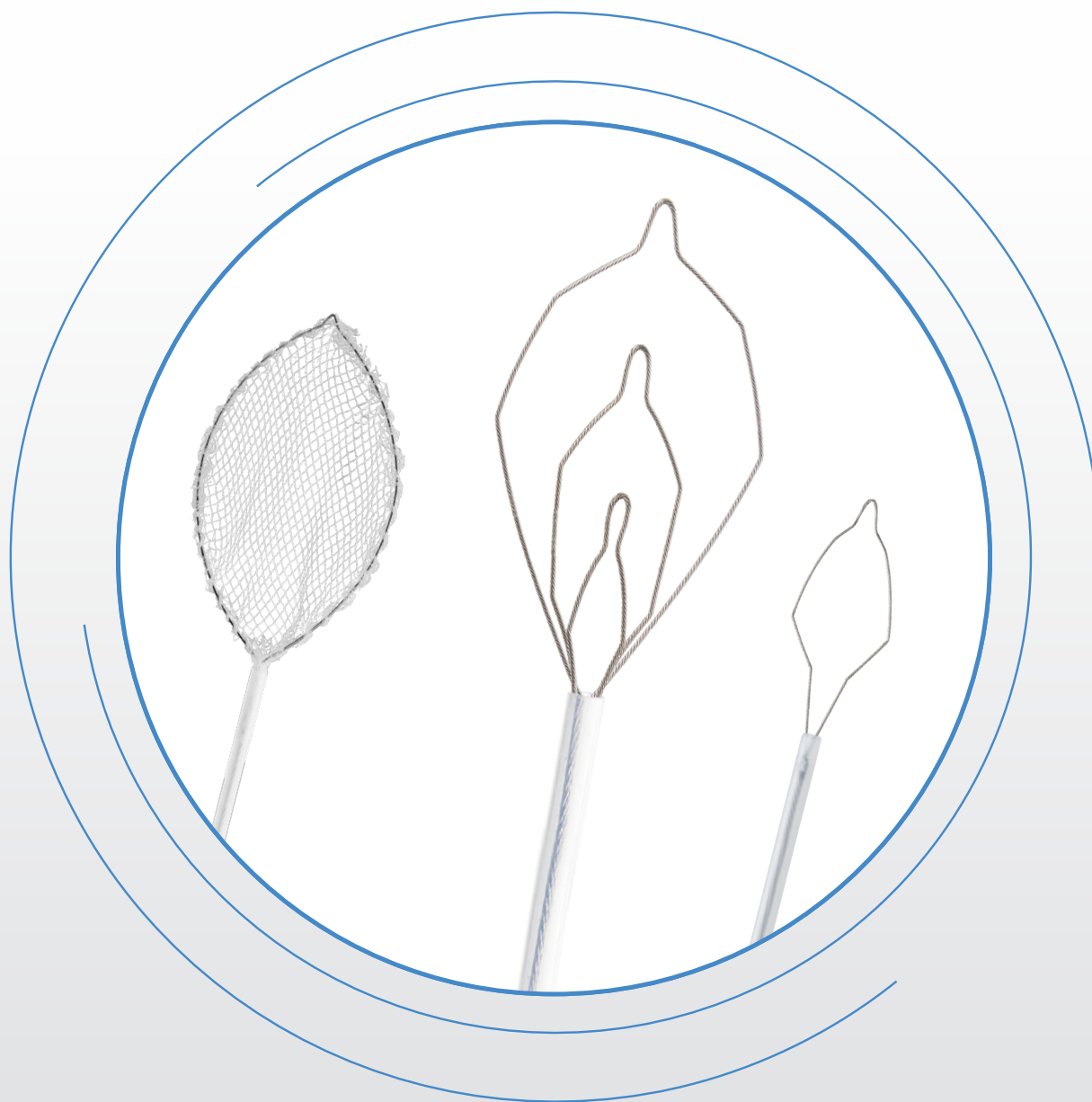

POLYPECTOMY:

From a dynamic hot and cold snare with a coined distal tip to our clinically proven cold snare, we offer a product solution to accommodate a variety of needs.

Polypectomy Solutions

TABLE OF CONTENTS:


Product Name	Page Number
SmartBand® Endoscopic Mucosal Resection Kit	2
Endoscopic Marker	3
COINTIP™ Snare	4
HISTOLOCK™ Resection Device	4
EXACTO™ Cold Snare	5
LARIAT™ Snare	6
ASCENDO™ Submucosal Lifting Agent	6
ROTH NET™ Retrievers	7
eTrap™ Polyp Trap	8
REVEAL™ Distal Attachment Cap	8
Product Information Chart	9

SmartBand®

Endoscopic Mucosal Resection Kit


SmartBand EMR Kit

- ▶ The SmartBand EMR Kit is indicated for endoscopic mucosal resection in the upper GI tract. This device is intended for single use only.
 - Consists of a ligation handle with universal connector, loading device, flush tube, Pentax adaptor, barrel with 5 latex free bands, and the SmartSnare™ EMR Hexagonal Snare.

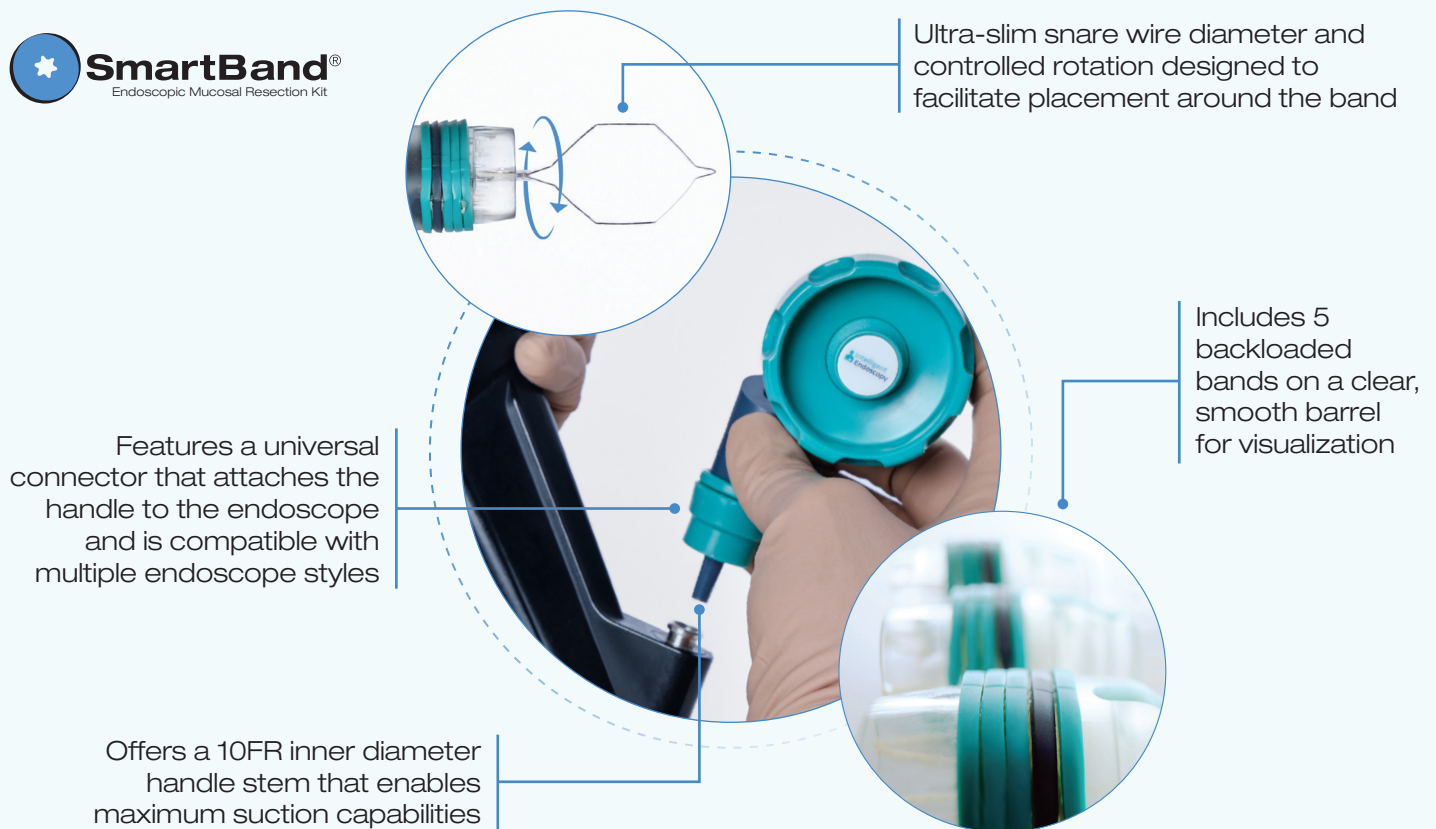

Key Features & Benefits...

- Offers a wide 360-degree visualization without the obstruction of bands
- Rotatable hex snare included
- Available in three barrel size options

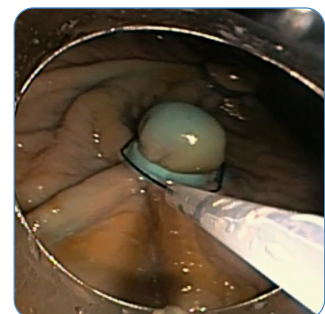

View of Smartband® EMR Kit procedure in porcine model

Endoscopic Marker

Endoscopic Marker

- ▶ The Endoscopic Marker is indicated for endoscopically marking lesions in the GI tract when the endoscopist anticipates the lesion will require surgical removal within 30 days.
- ▶ Provides a method of marking target tissue and lesions on the wall of a patient's stomach or colon during an endoscopic procedure. The mark provides a reference point for subsequent removal of the tissue.

Key Features & Benefits...

- *No shake necessary - Allows for stable suspension of the ink.*
- *Stable emulsion - Does not migrate resulting in dark and permanent marks.*
- *Each syringe is prefilled - Preparation is simplified with ready-to-use product. Each syringe comes with 5ml of the marking fluid.*
- *Sterile - The product is single use and sterile.*




COINTIP™ Snare / HISTOLOCK™ Resection Device


COINTIP Snare

► The COINTIP Snare features a dynamic braided wire configuration with a coined distal tip supporting hot and cold cutting capabilities specifically for complex polypectomies like endoscopic mucosal resection.

- *Distal tip flattened to 0.010" which is achieved by the manufacturing process of "coining"*
- *0.016" snare loop wire supports grip of tissue to capture target*
- *Catheter includes an internal strain relief for added support during resection to reduce buckling during tissue resection*


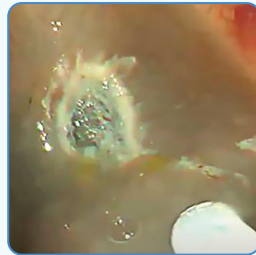


Resection sites after use of the COINTIP™ Snare.

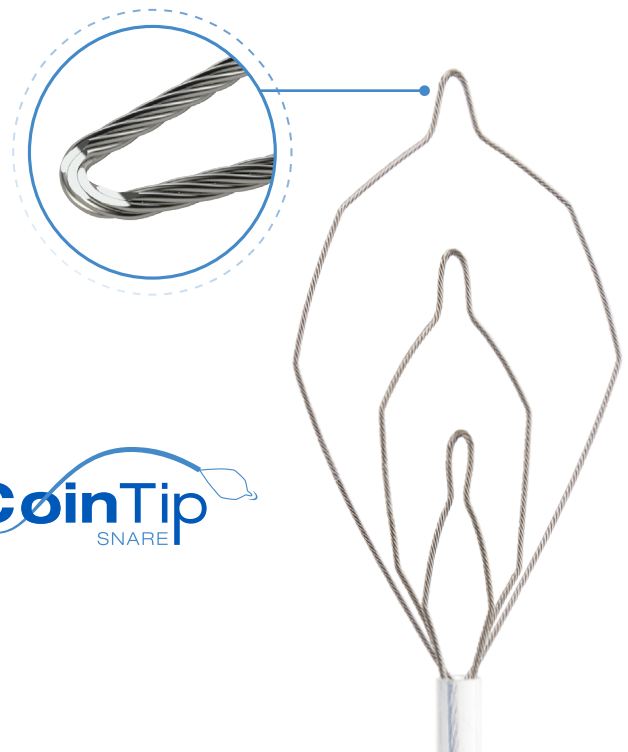


HISTOLOCK Resection Device

► The HISTOLOCK Resection Device features a stiff monofilament wire to help recruit and resect flat tissue during complex polypectomy procedures.

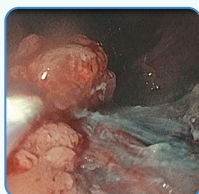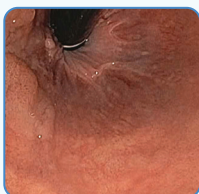

The HISTOLOCK™ Resection Device effectively resects a flat 4cm rectal polyp with precision.

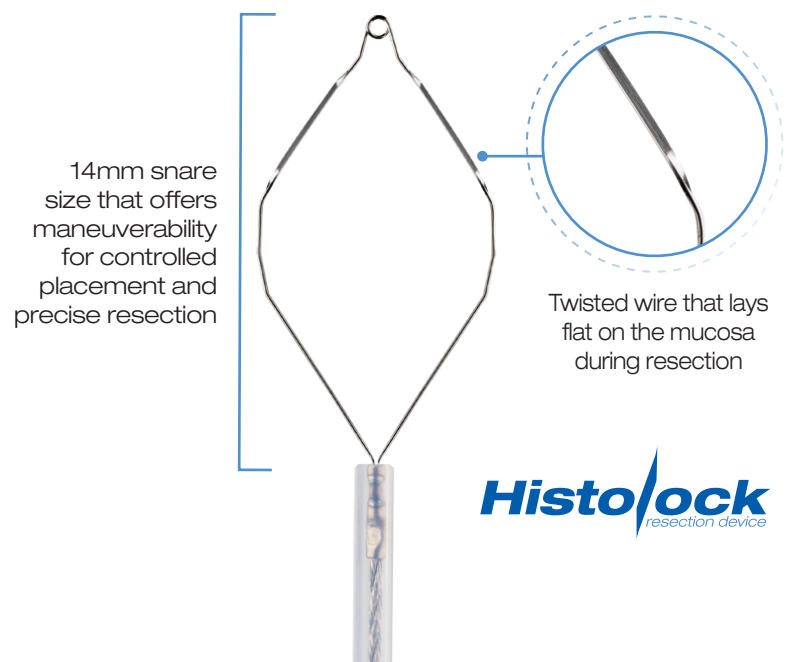

EXACTO™ Cold Snare

EXACTO Cold Snare


- ▶ Clinically proven to achieve a significantly higher rate of complete resection, the design of the EXACTO Cold Snare offers control and placement for a precise, clean cut when compared to traditional braided snares.^{1,2}
 - Enables a clean cut with a thin wire design (0.012"/0.30mm) while maintaining its shape and stiffness
 - Features a shield shape design that maximizes its width for control and placement
 - 9mm mini snare size that facilitates polyp tissue resection without the use of cautery
 - Offers a thin, stiff catheter that decreases buckling while aiding in placement



40mm laterally spreading tumor resected with the EXACTO™ Cold Snare



Clean, en bloc resection site of a subcentimeter polyp

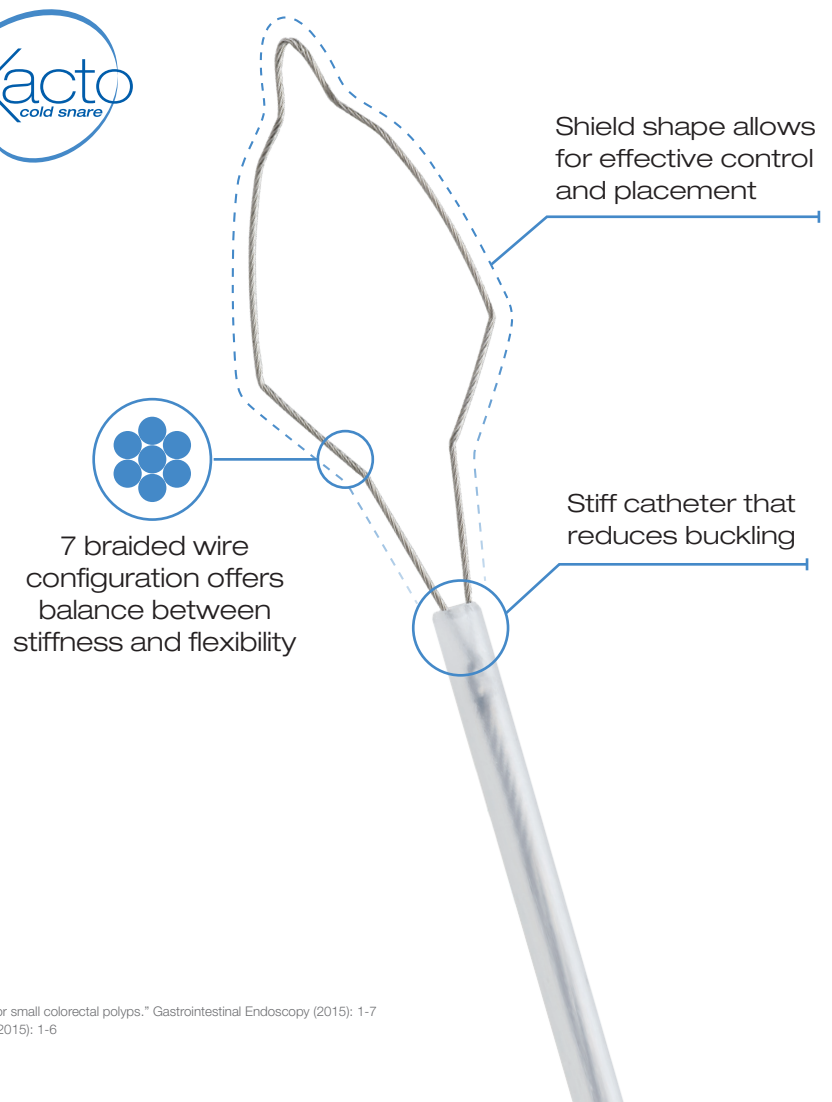

1. Horiuchi A, Hosoi K. "Prospective, randomized comparison of 2 methods of cold snare polypectomy for small colorectal polyps." *Gastrointestinal Endoscopy* (2015): 1-7
2. Din S, Ball A. "Cold snare polypectomy: Does snare type influence outcomes?" *Digestive Endoscopy* (2015): 1-6

LARIAT™ Snare /


ASCENDO™ Submucosal Lifting Agent

LARIAT Snare

▶ The multi-functional design of the LARIAT Snare helps address a variety of resection needs including piecemeal resection and multiple polyps.



- *Three distinct shapes with multiple sizes to reduce the need to exchange snares*
- *Snare loop design maintains its shape and integrity after multiple resections*

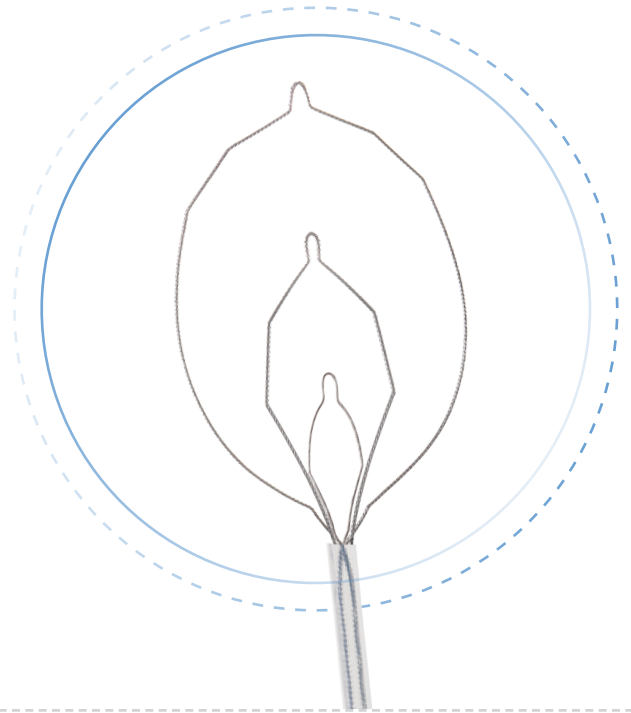

ASCENDO Submucosal Lifting Agent

▶ The ASCENDO Submucosal Lifting Agent helps to clarify the area for visualizing the margins of the target lesion and performing the resection procedure. It is designed to provide a submucosal cushion of height and duration.

- *Prefilled Syringe - Preparation is simplified with ready-to-use product*
- *Biocompatible - Consists of sodium hyaluronate, commonly referred to as hyaluronic acid, a biopolymer naturally present in the human body*




ROTH NET™ Retrievers


ROTH NET Retrievers

► From capturing large polyps to securing multiple polyp fragments, the ROTH NET Retriever is a device you can count on for consistent performance. Our all-encompassing line of retrieval nets...


- *Preserve polyp fragments and specimens for histopathology*
- *Allow for durability and continued visibility through the net*
- *Opens even in retroflexed position*

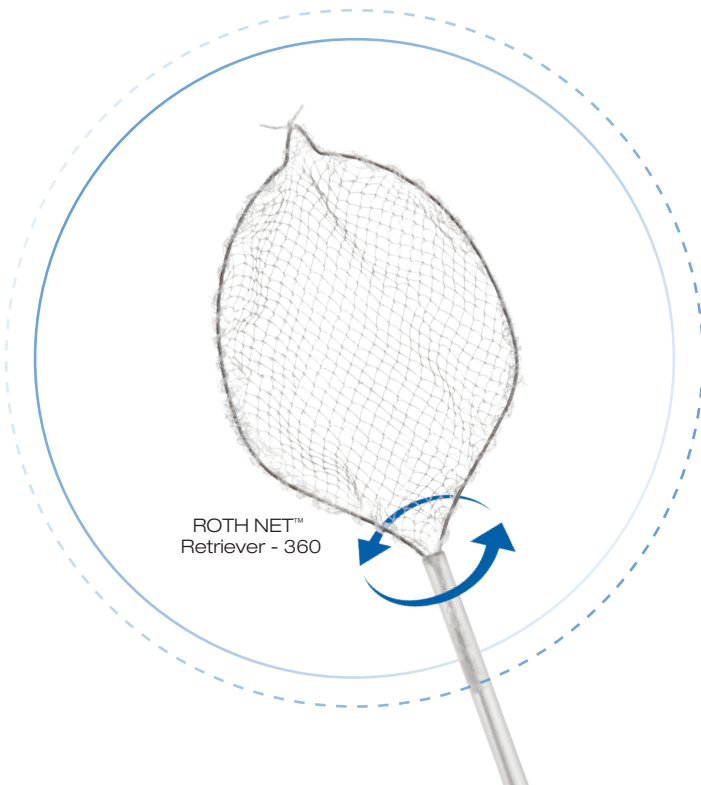

ROTH NET™
Retriever - 360

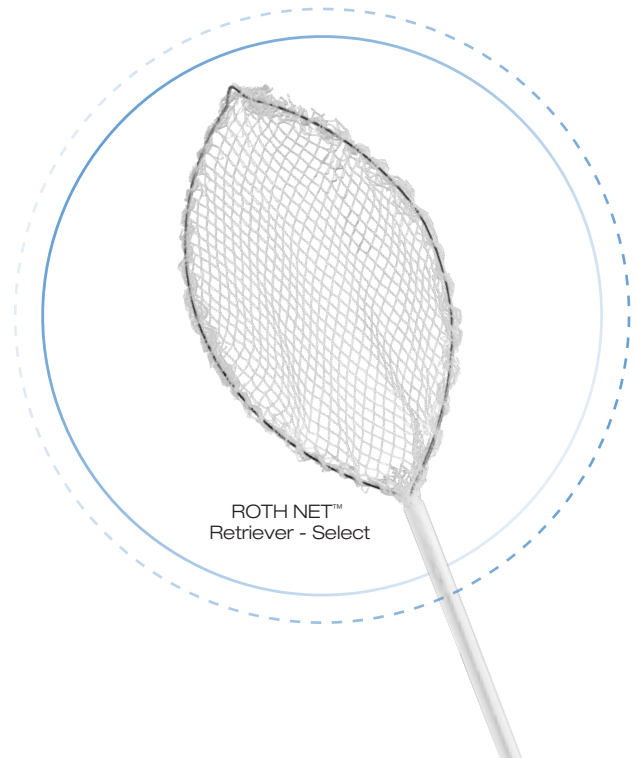


ROTH NET™
Retriever - Select

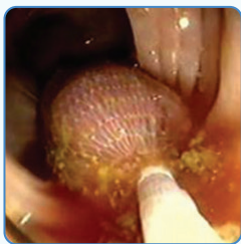

Retrieval of large polyp
with ROTH NET™ Retriever



Retrieval of multiple polyps/
polyp fragments with
ROTH NET™ Retriever

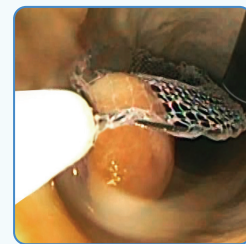

ROTH NET™ Retriever - select
offering durability and strength
across multiple passes

eTrap[®] Polyp Trap / REVEAL[®] Distal Attachment Cap

eTrap Polyp Trap

► Allows for retrieval of multiple specimens, while safeguarding clinicians from unnecessary exposure to biomaterials.


- *Collect multiple polyps with uninterrupted suction*
- *Simple transport and preservation of specimens*
- *Magnifying window allows direct visualization of the collected specimen*



REVEAL Attachment Cap

► The REVEAL Distal Attachment Cap is designed to facilitate the endoscope's view during polypectomy procedures.

- *Helps avoid 'red out' by maintaining a constant distance between the endoscope and mucosa*
- *Can be used with other polypectomy devices, including snares and needles*

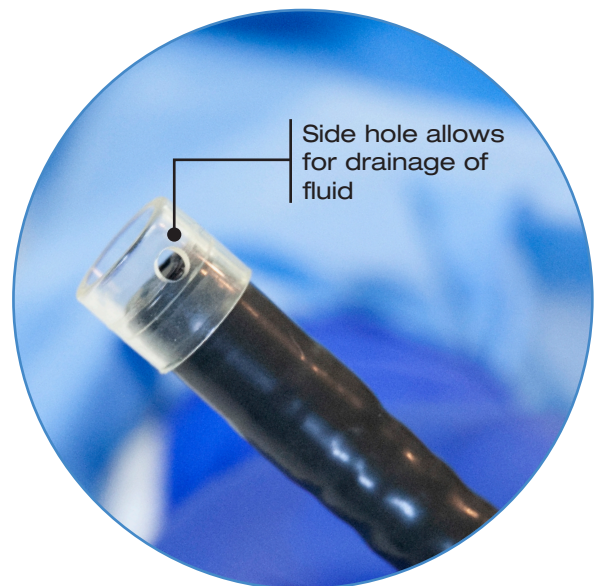

Product Information Chart

Product Number	Description	Barrel Size (mm)	Endoscopic Diameter Compatibility (mm)			Units/Box
SBEMRK94	SmartBand® EMR Kit 9.4mm Barrel Tip Diameter	9.4	8.6 – 11.6			EA
SBEMRK	SmartBand® EMR Kit 10.75mm Barrel Tip Diameter	10.75	8.6 – 11.6			EA
SBEMRK12	SmartBand® EMR Kit 12.0mm Barrel Tip Diameter	12.0	8.6 – 11.6			EA
SBEMRP94	SmartBand® EMR Pack 9.4mm Barrel Tip Diameter	9.4	8.6 – 11.6			EA
SBEMRP	SmartBand® EMR Pack 10.75mm Barrel Tip Diameter	10.75	8.6 – 11.6			EA
SBEMRP12	SmartBand® EMR Pack 12.0mm Barrel Tip Diameter	12.0	8.6 – 11.6			EA
Product Number	Description	Volume Per Syringe	Sterile			Units/Box
BX00711836	Endoscopic Marker	5ml	yes			10
BX00712934	ASCENDO™ Submucosal Lifting Agent	5ml	yes			10
Product Number	Description	Needle Gauge/Projection	Sheath Diameter (mm)	Length (cm)	Snare Size (mm)	Units/Box
BX00711018	SHORT THROW™ Snare - Mini	-	2.4	230	12.5 x 30	10
BX00711019	SHORT THROW™ Snare - Standard Oval	-	2.4	230	25 x 54	10
BX00711101	SHORT THROW™ Snare - Standard Oval	-	2.4	350	25 x 54	5
BX00711113	TRAXTION™ Hexagonal Snare	-	2.4	230	25 x 40	10
BX00711115	EXACTO™ Cold Snare	-	2.4	230	9 x 19	10
BX00711094	EXACTO™ Cold Snare	-	2.4	230	9 x 19	20
BX00711095	EXACTO™ Cold Snare	-	2.4	230	9 x 19	50
BX00711110	COINTIP™ Snare	-	2.5	230	Deltoid: 20 x 32 Hex: 10 x 27 Diamond: 6 x 19	10
BX00711117	HISTOLOCK™ Resection Device	-	2.4	230	14 x 27	10
BX00711871*	HISTOLOCK™ Resection Device	-	2.4	230	14 x 27	3
BX00711119	LARIAT™ Snare	-	2.4	230	Oval: 30 x 55 Hex: 10 x 28 Diamond: 6 x 15	10
BX00711096	LARIAT™ Snare	-	2.4	230		20
BX00711097	LARIAT™ Snare	-	2.4	230		50
BX00711120	ROTATOR™ Snare - Standard	-	2.4	230	25 x 54	10
BX00711121	ROTATOR™ Snare - Mini	-	2.4	230	12.5 x 30	10
Product Number	Description	Specimen Strainer Trays				Units/Box
BX00711099	eTrap™ Polyp Trap Visit steris.com/healthcare to view our full line of polyp traps	50 (2 per unit)				25
Product Number	Description	Endoscope Tip O.D.	Cap O.D.	Length	Packaging Color	Units/Box
BX00711770	REVEAL™ Distal Attachment Cap (Sterile)	8.9-9.5mm	11.35mm	4mm	Pink	10
BX00711771	REVEAL™ Distal Attachment Cap (Sterile)	9.5-9.9mm	11.8mm	4mm	Blue	10
BX00711772	REVEAL™ Distal Attachment Cap (Sterile)	9.8-10.5mm	12.4mm	4mm	Orange	10
BX00711773	REVEAL™ Distal Attachment Cap (Sterile)	10.9-11.8mm	13.4mm	4mm	Light Green	10
BX00711774	REVEAL™ Distal Attachment Cap (Sterile)	12.6-13.2mm	15.0mm	4mm	Yellow	10
BX00711775	REVEAL™ Distal Attachment Cap (Sterile)	13.4-13.9mm	15.7mm	4mm	Dark Green	10
Product Number	Description	Sheath Diameter (mm)	Length (cm)	Net Size (cm)	Qty/Box	
Tissue Retrieval Nets						
BX00711050	ROTH NET® Retriever - Standard	2.5	230	3 x 6	5	
BX00711151	ROTH NET® Retriever - Polyp	2.5	230	3 x 6	10	
BX00711162	ROTH NET® Retriever - Standard (Sterile)	2.5	230	3 x 6	10	
BX00715050	ROTH NET® Platinum® Retriever - Polyp	2.5	230	4 x 5.5	5	
BX00715060	ROTH NET® Platinum® Retriever - Polyp (Sterile)	2.5	230	4 x 5.5	5	
Foreign Body Retrieval Nets						
BX00711150	ROTH NET™ Retriever - Standard	2.5	230	3 x 6	5	
BX00711052	ROTH NET™ Retriever - Standard	2.5	230	3 x 6	10	
BX00711053	ROTH NET™ Retriever - Standard	2.5	160	3 x 6	5	
BX00711169	ROTH NET™ Retriever - Select (Sterile)	2.5	230	3 x 6	10	
BX00711198	ROTH NET™ Retriever - 360 (Sterile)	2.5	230	3 x 5	10	
BX00711200	ROTH NET™ Retriever - 360	2.5	230	3 x 5	10	
BX00711155	ROTH NET™ Platinum Retriever - Food Bolus	2.5	160	4 x 5.5	5	
BX00711165	ROTH NET™ Retriever - Maxi (Sterile)	3.0	230	3 x 6	5	
BX00715060	ROTH NET™ Platinum® Retriever - Universal (Sterile)	2.5	230	4 x 5.5	5	
BX00715051	ROTH NET™ Platinum® Retriever - Universal	2.5	230	4 x 5.5	10	
Specialty Retrieval Nets						
BX00711152	ROTH NET™ Retriever - Endoscope	2.5	350	3 x 6	5	
BX00711158	ROTH NET™ Retriever - Mini Polyp	2.8	230	2.5 x 4.5	10	
BX00711159	ROTH NET™ Retriever - Mini Polyp (Sterile)	2.8	230	2.5 x 4.5	5	
BX00711168	ROTH NET™ Retriever - Mini (Sterile)	1.8	160	2 x 4.5	5	

*Not for sale in the U.S.



ROTH NET® and ROTH NET® Platinum® is a registered trademark of STERIS Endoscopy.

©2023 STERIS Corporation. All rights reserved.
All company and product names are trademarks of STERIS Corporation,
its affiliates or related companies, unless otherwise noted.

Follow us:    



5960 Heisley Road
Mentor, OH 44060-1834 – USA

Customer Service: 800-548-4873

Fax: 440-639-4495

www.steris.com