

SmartBand® Ligation System

Scientifically designed to improve ligation performance



Product Specifications


Used to endoscopically ligate esophageal varices at or above the gastroesophageal junction and to ligate internal hemorrhoids.




High Retention Bands (HRB) deliver the band performance needed to treat the patient with a 24-month shelf life.



Predictable band deployment using modern, ergonomic science with just a natural twist of the wrist to deploy each band and prepare the next band for deployment.

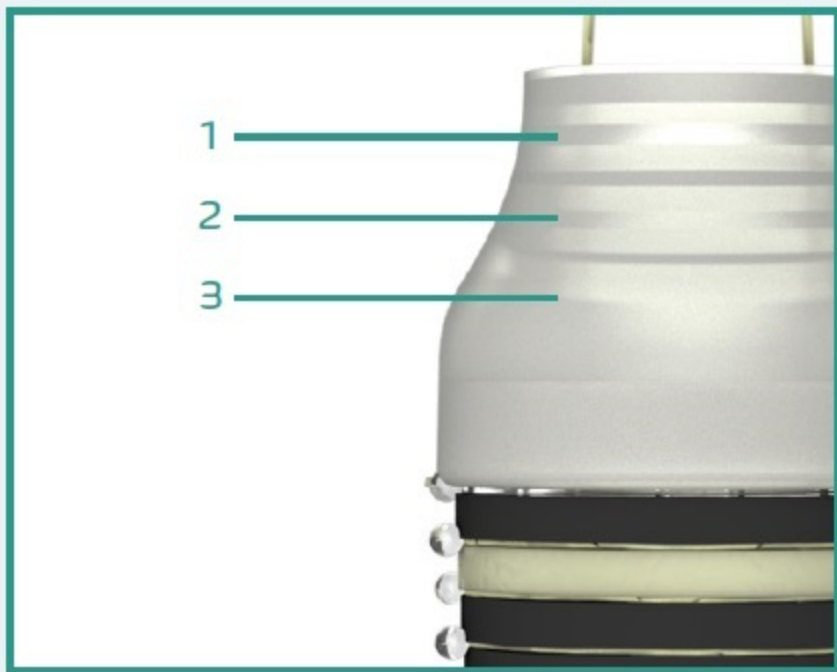

3-Intelligent Ring technology creates a customized and secure endoscope and barrel connection during intubation and extubation.

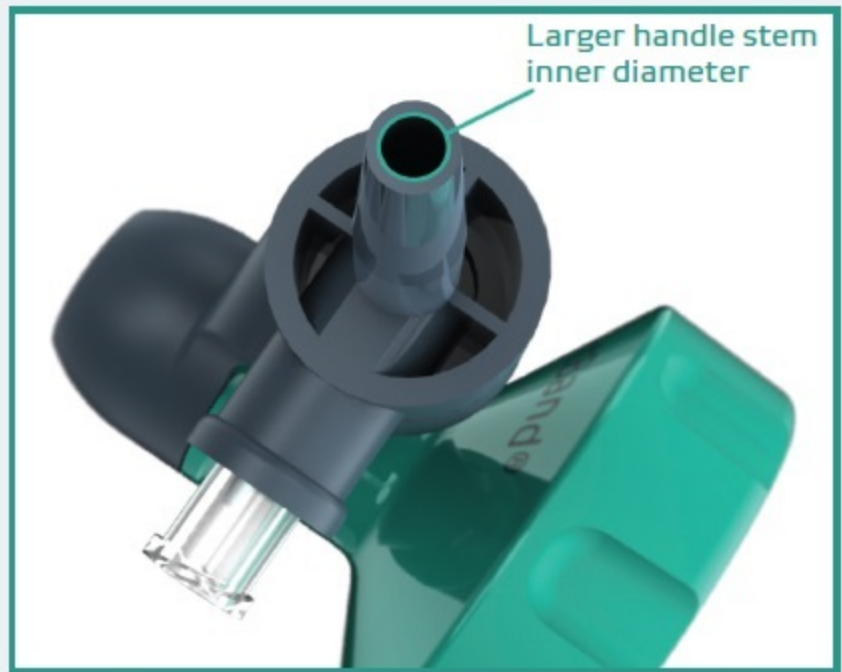

The industry's largest handle inner stem diameter of 9.8FR ensures maximum clearance for suction and flushing capabilities throughout the procedure.



A dedicated flush port and connector tube makes flushing the endoscope channel fast and easy keeping it clean and free of debris.



The SmartBand Ligation Pack, "SmartPack", is now available when more bands are needed to complete the procedure.

SmartBand Ligation Kit

Components: Ligation Handle, Loading Device, Connector Tube and Deployment Cord with Barrel of Bands

| Order Number | Reference Part Number | Barrel Diameter (mm) | Number of Bands | Band Material | Shelf Life (months) |
|--------------|-----------------------|----------------------|-----------------|---------------|---------------------|
| IE11002 | SLK-6 | 8.6 - 11.6 | 6 | Latex | 24 |

This device is supplied non-sterile and is disposable, single-use only

SmartBand Ligation Pack

Components: Deployment Cord with Barrel of Bands

| Order Number | Reference Part Number | Barrel Diameter (mm) | Number of Bands | Band Material | Shelf Life (months) |
|--------------|-----------------------|----------------------|-----------------|---------------|---------------------|
| IE11001 | SLP-6 | 8.6 - 11.6 | 12 | Latex | 24 |

This device is supplied non-sterile and is disposable, single-use only

Customer Service:

orders@intelligentendoscopy.com

Phone: +1 (336) 608-4375
Toll Free: +1 (844) 344 9376
Fax: +1 (336) 893-7268

www.intelligentendoscopy.com



4740 Commercial Park Court
Suite 1
Clemmons, NC 27012

